

PRECISION SPINE REFORM[®] PEDICLE SCREW SYSTEM

STRENGTH, STABILITY AND EFFICIENCY

Reform[®] is a comprehensive pedicle screw system that is designed to meet the varying requirements of degenerative, trauma and deformity procedures. Reform features a cobalt chrome tulip, a titanium triple lead proximaly tapered thread and titanium and cobalt chrome rods to deliver strength, stability and efficiency to all thoracolumbar constructs. Reduction and uniplanar screw options, along with a full line of hooks, dominoes and offsets, complete the system to simplify the procedure and accommodate individual patient anatomy.



COBALT CHROME TULIP



UNIPLANAR REDUCTION



TRIPLE LEAD THREAD



Discover the Difference



REFORM® PEDICLE SCREW SYSTEM



Cobalt Chrome Tulip
reduces risk of head splay
with DVR maneuvers

**Square Threaded
Locking Cap**
geometry reduces risk
of cross threading



Hexalobular Drive Fittings
on Bone Screw and Locking
Cap reduces incidence of
toggle and stripping




Polyaxial Screws
provide 30° angulation in all
planes for intraoperative flexibility

Proximal Tapered Thread
design increases bone/screw
interface enhancing pull-out
strength



**Reduction Tower
Quick Connect**
allows for slip-on
insertion over
the tulip

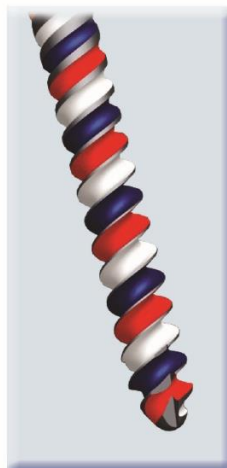

Triple Lead Thread
allows for efficient
screw delivery

REFORM® SIZING

Screw Sizes

4.5	25-45mm (5mm)
5.5	30-50mm (5mm)
6.5	30-55mm (5mm)
7.5	35-55mm (5mm)
8.5	40-60mm (5mm)
	70 & 80mm
9.5	60-80mm (10mm)

Rod Diameter

5.5mm

Ti Straight

80, 100, 120, 200, 300, 400
& 500mm

Ti Curved

35-80mm (5mm)
90-120mm (10mm)

CoCr Straight

80, 100, 120, 200, 300, 400
& 500mm

CoCr Curved

35-80mm (5mm)
90-120mm (10mm)

Precision Spine, Inc.
2050 Executive Drive, Pearl, MS 39208
Customer Service: 1.888.241.4773
Phone: 601.420.4244
Toll Free: 877.780.4370
Fax: 601.420.5501
www.precisionspineinc.com



CE0086